

# Infinity<sup>®</sup> LX

Medical Image Server



**CODONICS**

A new dimension in image storage

# Infinity<sup>®</sup> LX

## Medical Image Server

### Overview

Infinity LX Medical Image Server is a DICOM storage appliance that provides a simple solution for storing and retrieving medical images. Applications such as direct modality storage of large DICOM data sets or storage of digital images in an environment without a centralized image server are ideal for Infinity. Infinity makes the power of complicated, rack-mounted servers available to users without special IT skills and is designed for expansion. With an intuitive interface and many pre-configured features, users can set up and maintain their own Infinity image servers. Infinity is backed by Codonics global service and support network to provide years of worry-free operation. When simple, affordable image storage is needed, Infinity is the answer.

Infinity's unique design delivers performance and reliability in an easy to use package. This makes Infinity an ideal storage solution for large, temporary datasets from CT, MR and other DICOM modalities. Thin slices, raw data, 3D reconstructions and multi-frame studies can be quickly saved and retrieved for future reconstruction or review. Automatic image management removes older studies tagged for temporary storage when space is needed. Infinity provides radiology departments with a simple solution that saves time and money when storing data not requiring long term archive on centralized servers.

### Ease of Setup and Use

Infinity includes many preconfigured settings for common applications so users can begin utilizing the system quickly. While Infinity's default configuration can be up and running within minutes, it is also highly customizable with a rich set of features including DICOM store, query/retrieve, storage commit, on-the-fly image compression, workstation access control and temporary or permanent storage of images. Infinity provides a comprehensive set of system notifications and warnings using web, email, RSS and visual/audible alerts built into the system. You take your data seriously, so does Infinity.

### Features / Benefits

- Compact, all-in-one design provides an affordable solution for DICOM image storage
- High performance architecture for fast DICOM file transfers
- Optimized for large DICOM data sets including DICOM files over 2GB
- Highly redundant hardware design including RAID, power supplies and cooling fans
- Simple system setup with pre-configured settings for non-technical personnel
- Optional Basic Web Viewer for quick JPEG viewing and optional Advanced Remote viewer for detailed DICOM images. These viewers feature multiple user access with easy deployment and maintenance.
- Optional JPEG lossless image compression increases storage capacity up to 2.5x
- Configurable access control with audit trails to assist with HIPAA compliance

### Specifications



User Interface:	Web and Windows <sup>®</sup> remote desktop
Protocols:	DICOM 3.0 (store, storage commit, query/retrieve)
Storage:	7.7 TB usable storage (3.5" SATA-2 hard drives)
Performance:	50 Mbytes/second (DICOM store)
Redundancy:	RAID 1 (mirroring) for operating software and database RAID 5 for image data Hot swap hard drives Hot swap power supply (1+1) Hot swap internal cooling fans Hot spare online hard drive with copy-back (optional)
Compression:	JPEG lossless or lossy compression of DICOM images (optional)
Alarms:	Drive failure, cooling fan failure, power supply failure, over-temperature warning, storage full, UPS status
Notifications:	Web, Email, RSS, built-in visual and audible alarms
Support:	Remote diagnostic/support via Internet (LogMeln Rescue <sup>®</sup> )
UPS:	Fully integrated UPS support with auto-shutdown
Power:	Universal Input: 100-240 VAC, 50/60 Hz, 650 Watts
Dimensions:	8.7" (22 cm) W, 26" (66 cm) D, 16.9" (43 cm) H
Weight:	92 lbs. (42 kg.) with drives installed



### Web System Status Interface



### Advanced Remote Viewer

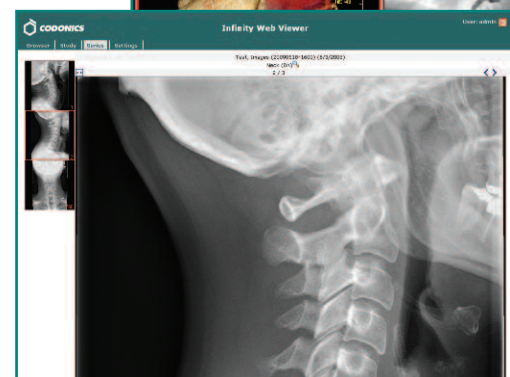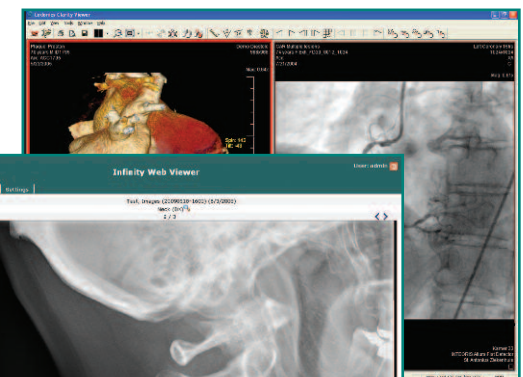

### Basic Web Viewer



Call Codonics today at **800-444-1198** or visit [www.codonics.com](http://www.codonics.com) for more information.

All registered and unregistered trademarks are the property of their respective owners. Specifications subject to change without notice. Printed in the U.S.A. Copyright © 2008-2009 Codonics, Inc. 9/2009

17991 Englewood Drive  
Middleburg Heights, OH 44130 USA  
(440) 243-1198  
(440) 243-1334 Fax  
Email [info@codonics.com](mailto:info@codonics.com)  
[www.codonics.com](http://www.codonics.com)

900-442-004.04